


PERKINS HAND HELD APPLANATION TONOMETER MK 2 PLUS

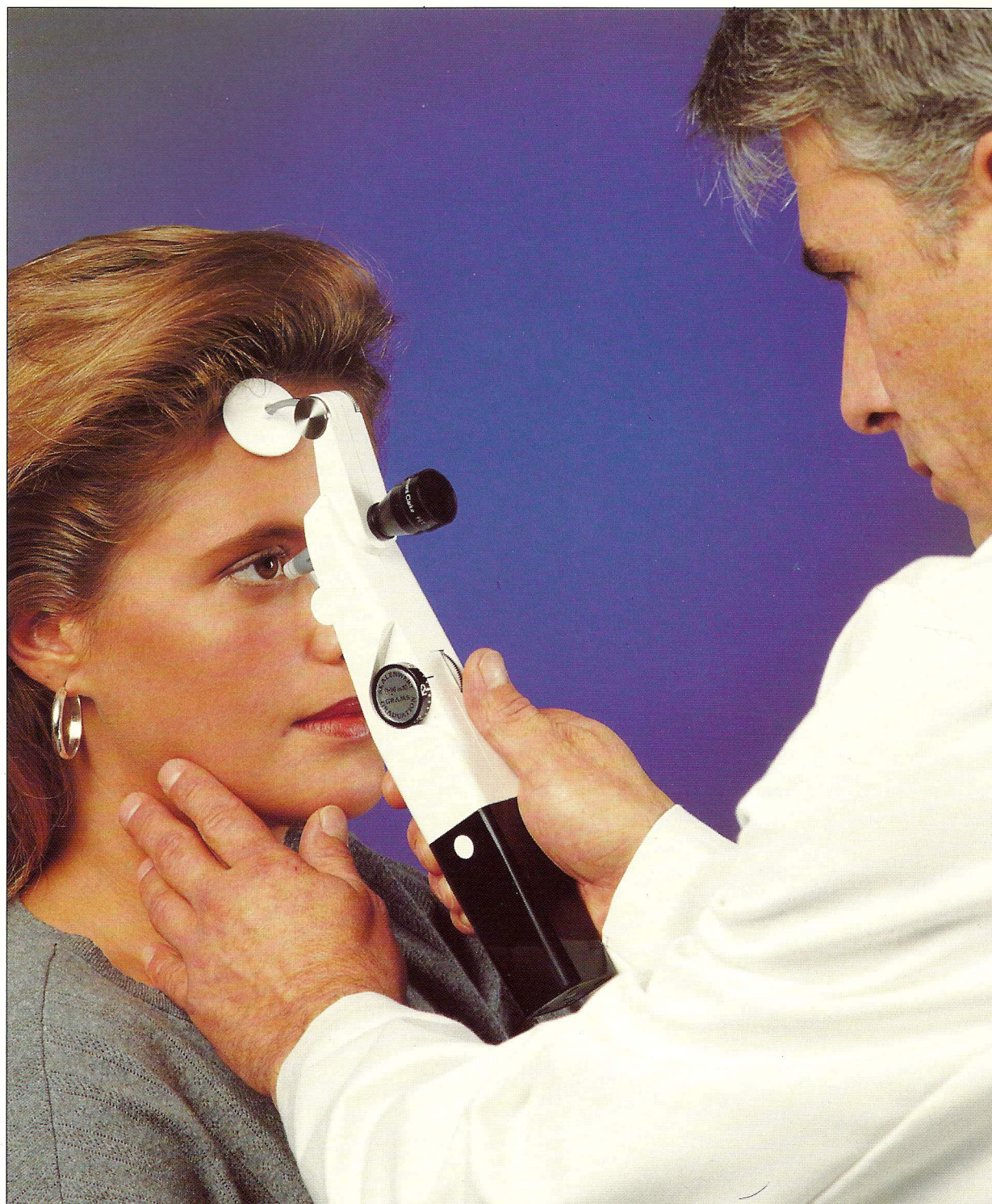

Clement Clarke



A HAAG-STREIT COMPANY

THE PERKINS HAND HELD APPL



The Perkins Tonometer has been established for over 20 years as a "world standard" in hand-held applanation tonometers.

Based on Goldmann principles and incorporating the original prisms from Haag-Streit, the Perkins Tonometer now has the advantage of combining well established benefits with new features to provide an instrument with even better all round performance.

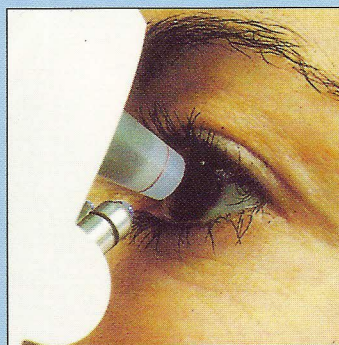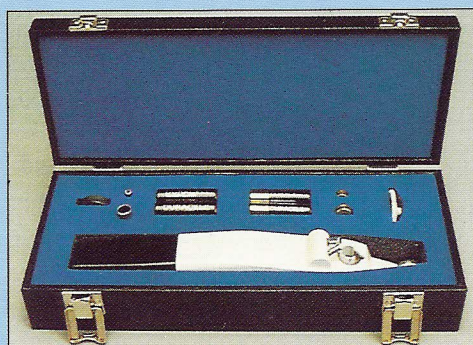

- **Primary Choice**

- The Perkins Tonometer can be used where other types of tonometer cannot! Post-operative and bed-ridden patients in particular are difficult to examine by traditional or non-contact methods and so the Perkins really is a primary choice of instrumentation for many people.

- **Goldmann Standards**

- The Perkins Tonometer uses the same prisms and test standards as the Goldmann Tonometer by Haag-Streit A.G.

- **Portability**

- The Perkins Tonometer is supplied in a compact case which is ideal for domiciliary visits.

- **Improved illumination**

- Twin lamps provide superior fluorescence which greatly enhance the familiar "semi-circles" and make operation of the instrument that much easier.

- **Versatility**

- Precision counter-balanced movements allow readings to be taken with the patient in any position and without affecting accuracy.


- **Accuracy**

- A re-designed operating mechanism increases sensitivity and every instrument is supplied with a signed and dated Certificate of Accuracy.

ANATION TONOMETER Mk 2 PLUS

- **Easy Maintenance**

- Spares and accessories are supplied as standard and the battery handle detaches easily for quick installation or removal of batteries. Four commercial AA size batteries may be used. Polarity is clearly marked on the bottom of the battery handle.

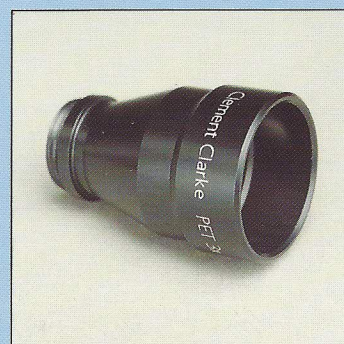


The Perkins Examination Telescope (P.E.T.) is a NEW accessory which allows MAGNIFIED viewing of the fluorescein "semi-circles" at arms length - making examinations easier for both you and the patient.

It also allows two people to view simultaneously and is an excellent teaching aid.

If you already use a Perkins Mk2 Tonometer the P.E.T. is fully compatible and mounts easily into the existing eyepiece aperture.

** New Release*

**PERKINS EXAMINATION
TELESCOPE (P.E.T.)**

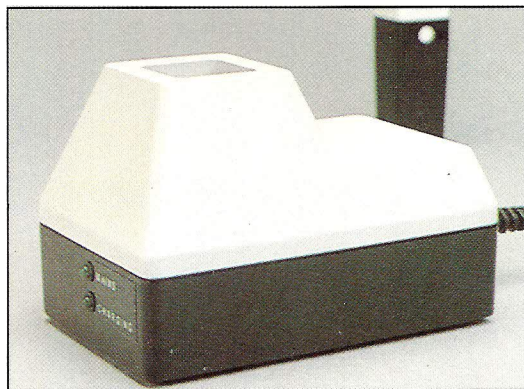

- **Rechargeable Version**

- A rechargeable battery handle is interchangeable with the standard handle and can be left to charge overnight in a custom- designed recharging unit.

ORDERING DETAILS

Part No Description

5805001 Perkins Hand Held Applanation Tonometer
Mk2 with carrying case and spares ☐

5805380 Recharging Unit ☐

5805490 Rechargeable Battery Handle ☐

5814004 Perkins Examination Telescope (P.E.T.) ☐

DISTRIBUTED BY

